
CE
0120

FORESEE LINE
S P R L

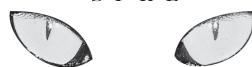RUE DU MANY, 111
4100 SERAING - BELGIUM
TEL. : 32(0)475 86 86 35
FAX : 32(0)4 388 15 45
e-mail : info@foreseeline.com

"STAR'TEX" CONDOMS
Instructions

DESCRIPTION AND INDICATIONS

Transparent lubricated condom with teat.
"STAR'TEX" condoms are made from natural latex.
Lubricant: Non-spermicidal dimethylsilicone.
Each condom is electronically tested and individually packed.

If used correctly, "STAR'TEX" condoms are an effective contraceptive and a safe way of combating the spread of sexually transmitted diseases.
However, no condom can guarantee complete protection.

INSTRUCTIONS

- Carefully tear the packet where indicated, taking care not to catch the condom on any sharp or pointed objects (nails, rings etc.)
- The condom should be placed on the erect penis before any contact with your partner.
- Hold the condom by the teat (rounded, sealed end) between your thumb and index finger in order to expel any air.
- Unroll the condom fully over the erect penis.
- After ejaculation and while the penis is still erect, withdraw the penis whilst holding the condom at the base.
Avoid all contact between the penis, the used condom and the body of your partner.
- After use, wrap the used condom in paper and throw it in the bin.
- Condoms should only be used once. Never use the same condom more than once.

GENERAL INFORMATION

LUBRICATION: if you require additional lubrication, use a water-based lubricant, preferably "STAR'TEX" gel. Never use vaseline, handcream etc.

COMPATIBILITY: if you are having medical treatment for your reproductive organs, consult your doctor or pharmacist to make sure that the medication is compatible with the condom.

ALLERGIES: if you are allergic to latex, consult a doctor.

STORAGE: store at room temperature, never leave in direct sunlight or very hot places.

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